




JDENTALCARE

just smile

JD ICON PLUS **T**

The power of T design

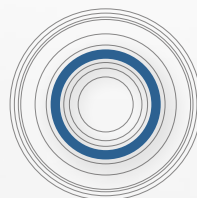


conical connection



Made in Italy


PRODUCT FEATURES

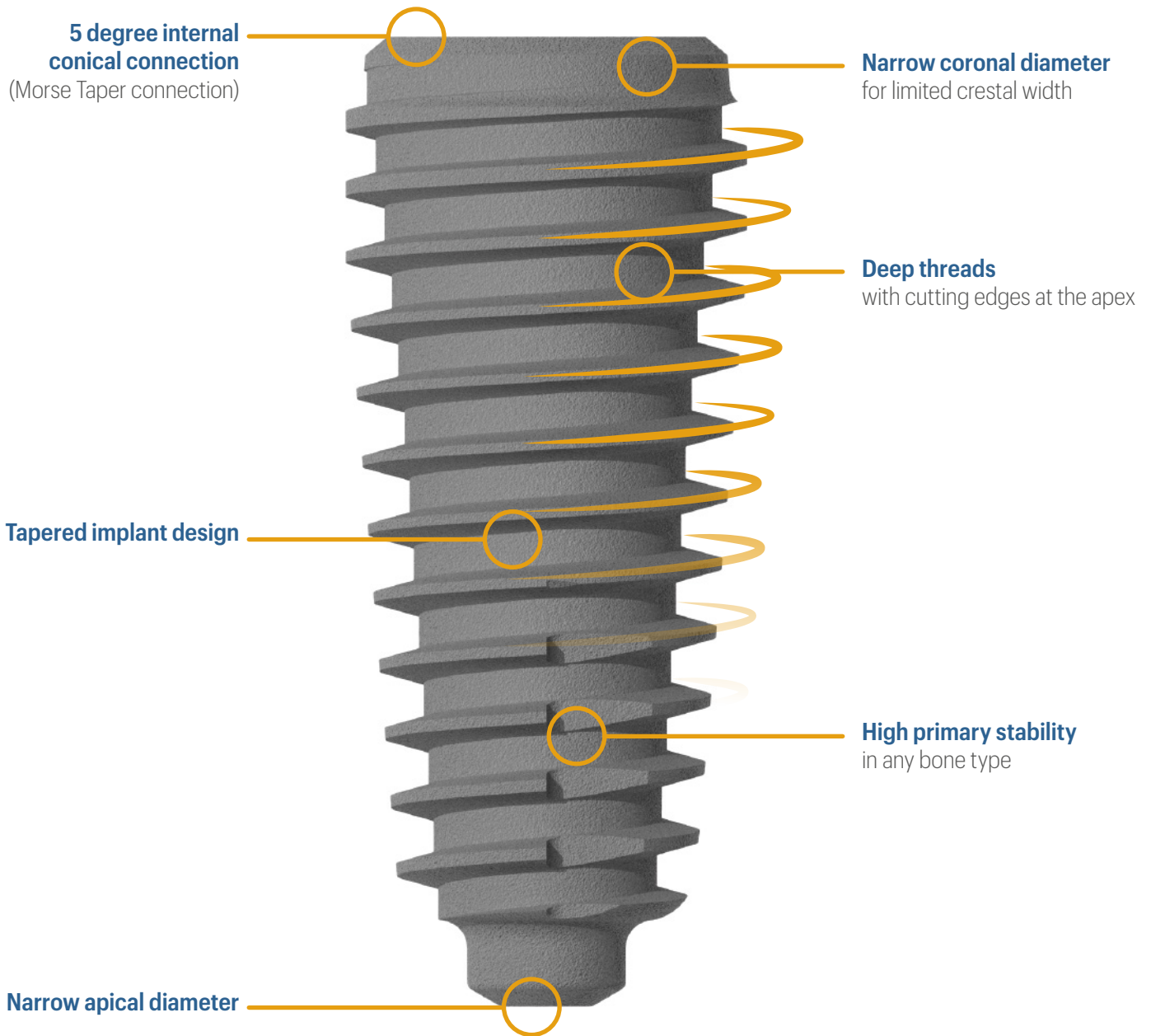

PRODUCT SPECIFICATIONS

The JDIcon Plus T dental implants are available in the diameters and lengths shown in the following table:

IMPLANT DIAMETER	LENGTH	TIP DIAMETER	PLATFORM DIAMETER	ABUTMENT INTERFACE
Ø 3,5	10, 11.5, 13	1.6	3.4	3.1
Ø 4,0	8.5	1.75	3.5	3.1
	10, 11.5, 13		3.6	
Ø 4,5	8.5, 10, 11.5, 13	1.75	3.6	3.1
Ø 5,0	8.5, 10, 11.5, 13	2.2	4.0	3.1
Ø 6,0	8.5, 10, 11.5	3	5.0	3.1
Ø 7,0	7, 8.5	2.5	5.0	3.1
Ø 8,0	7, 8.5	2.8	5.0	3.1

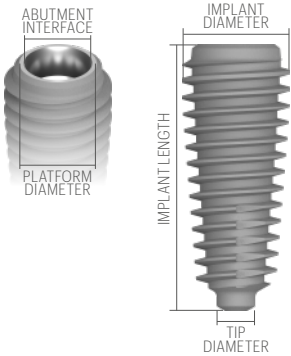

Note: all measurements in mm

PRODUCT CATALOGUE

Note: all measurements of diameters and lengths are in mm

Implants:


IM37100-0	JDIcon Plus Ø 3.5 L 10			
IM37115-0	JDIcon Plus Ø 3.5 L 11.5			
IM37130-0	JDIcon Plus Ø 3.5 L 13			
ICT40080.	JDIcon Plus Ø 4.0 L 8.5			
ICT40100.	JDIcon Plus Ø 4.0 L 10			
ICT40115.	JDIcon Plus Ø 4.0 L 11.5			
ICT40130.	JDIcon Plus Ø 4.0 L 13			
ICT45080.	JDIcon Plus Ø 4.5 L 8.5			
ICT45100.	JDIcon Plus Ø 4.5 L 10			
ICT45115.	JDIcon Plus Ø 4.5 L 11.5			
ICT45130.	JDIcon Plus Ø 4.5 L 13			
ICT50080.	JDIcon Plus Ø 5.0 L 8.5			
ICT50100.	JDIcon Plus Ø 5.0 L 10			
ICT50115.	JDIcon Plus Ø 5.0 L 11.5			
ICT50130.	JDIcon Plus Ø 5.0 L 13			
ICT55080.	JDIcon Plus Ø 6.0 L 8.5			
ICT55100.	JDIcon Plus Ø 6.0 L 10.0			
ICT55115.	JDIcon Plus Ø 6.0 L 11.5			
ICT70070.	JDIcon Plus Ø 7.0 L 7.0			
ICT70080.	JDIcon Plus Ø 7.0 L 8.5			
ICT80070.	JDIcon Plus Ø 8.0 L 7.0			
ICT80080.	JDIcon Plus Ø 8.0 L 8.5			

Note: All the prosthetic components are the same as JDIcon Plus implant line

DRILLING PROTOCOLS

It is recommended to adhere to the indications of the following drilling sequence to ensure optimal primary stability of the implants.

	HEALED BONE		POST EXTRACTIVE BONE	
IMPLANT DIAMETER	SOFT BONE	MEDIUM-DENSE BONE	SOFT BONE	MEDIUM-DENSE BONE
site preparation in maxilla				
Ø 3,5	2,0 2,4 2,8 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 at the entrance	2,0 2,4 2,8 3,2 up to the 1 st laser mark L 6mm
Ø 4,0	2,0 2,4 2,8 3,2 at the entrance	2,0 2,4 2,8 3,2 up to the 1 st laser mark L 6-8mm 3,6 at the entrance	2,0 2,4 2,8 3,2 at the entrance	2,0 2,4 2,8 3,2 at the entrance
Ø 4,5	2,0 2,4 2,8 3,2 3,6 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 3,6 up to the 1 st laser mark L 6mm 4,0 at the entrance	2,0 2,4 2,8 3,2 3,6 at the entrance	2,0 2,4 2,8 3,2 3,6 up to the 1 st laser mark L 6mm
Ø 5,0	2,0 2,4 2,8 3,2 3,6 4,0 at the entrance	2,0 2,4 2,8 3,2 3,6 4,0 up to the 1 st laser mark L 6mm 4,4 at the entrance 4,8 at the entrance	2,0 2,4 2,8 3,2 3,6 at the entrance	2,0 2,4 2,8 3,2 3,6 up to the 1 st laser mark L 6mm 4,0 at the entrance
Ø 6,0	2,0 2,4 2,8 3,2 3,6 4,0 4,4 up to the 1 st laser mark L 6mm 4,8 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 3,6 4,0 4,4 up to the 1 st laser mark L 6mm 4,8 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 3,6 4,0 4,4 up to the 1 st laser mark L 6mm 4,8 at the entrance	2,0 2,4 2,8 3,2 3,6 4,0 up to the 1 st laser mark L 6mm 4,4 at the entrance 4,8 at the entrance
site preparation in mandible				
Ø 3,5	2,0 2,4 2,8 3,2 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 up to the 2 nd laser mark L 8mm 3,6 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 at the entrance	2,0 2,4 2,8 3,2 up to the 2 nd laser mark L 8mm
Ø 4,0	2,0 2,4 2,8 3,2 3,6 up to the 2 nd laser mark L 8mm	2,0 2,4 2,8 3,2 3,6 4,0 up to the 2 nd laser mark L 8mm 4,4 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 at the entrance	2,0 2,4 2,8 3,2 3,6 at the entrance
Ø 4,5	2,0 2,4 2,8 3,2 3,6 at the entrance 4,0 at the entrance	2,0 2,4 2,8 3,2 3,6 4,0 4,4 4,8 at the entrance	2,0 2,4 2,8 3,2 3,6 at the entrance	2,0 2,4 2,8 3,2 3,6 4,0 at the entrance
Ø 5,0	2,0 2,4 2,8 3,2 3,6 at the entrance 4,0 at the entrance	2,0 2,4 2,8 3,2 3,6 4,0 4,4 4,8 at the entrance	2,0 2,4 2,8 3,2 3,6 at the entrance	2,0 2,4 2,8 3,2 3,6 4,0 at the entrance
Ø 6,0	2,0 2,4 2,8 3,2 3,6 4,0 4,4 4,8 up to the 1 st laser mark L 6mm 5,4 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 3,6 4,0 4,4 4,8 5,4 up to the 2 nd laser mark L 8mm	2,0 2,4 2,8 3,2 3,6 4,0 at the entrance 4,4 at the entrance 4,8 at the entrance	2,0 2,4 2,8 3,2 3,6 4,0 up to the 1 st laser mark L 6mm 4,4 up to the 1 st laser mark L 6mm 4,8 up to the 1 st laser mark L 6mm

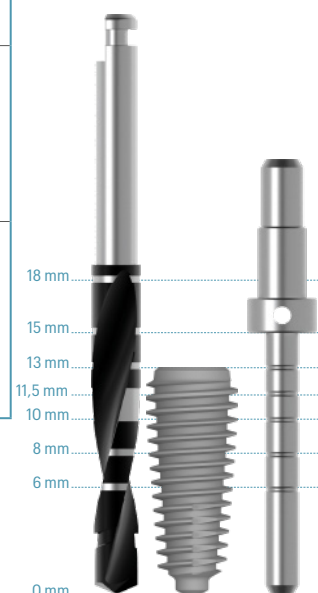


Note: all measurements in mm

DRILLING PROTOCOLS - Ø 7.0 mm / Ø 8.0 mm

Implant diameters 7 mm and 8 mm are recommended for posterior maxillary sites and in molar post extractive bone.

IMPLANT DIAMETER	SOFT BONE	MEDIUM-DENSE BONE
Ø 7,0	2,0	2,0
	2,4	2,4
	2,8	2,8
	3,2	3,2
	3,6	3,6
	4,0	4,0
	4,4	4,4
	4,8 up to the 1 st laser mark L 6mm	4,8 up to the 1 st laser mark L 6mm
Ø 8,0	2,0	2,0
	2,4	2,4
	2,8	2,8
	3,2	3,2
	3,6	3,6
	4,0	4,0
	4,4	4,4
	4,8 up to the 1 st laser mark L 6mm	4,8 up to the 1 st laser mark L 6mm
	5,4 up to the 1 st laser mark L 6mm	5,4 up to the 1 st laser mark L 6mm

Note: all measurements in mm

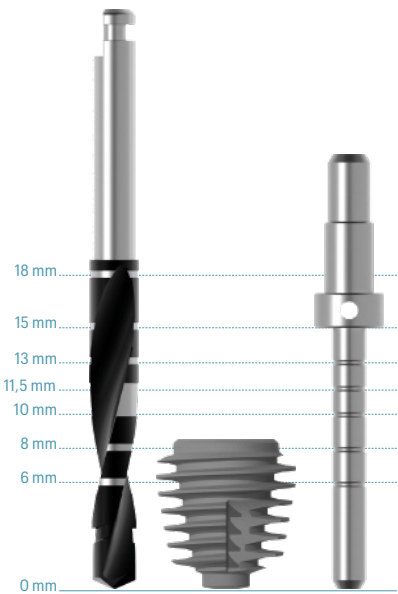



JDENTALCARE

www.jdentalcare.com